



## Otosure

# PC-based automatic screening audiometer

For fast, accurate hearing tests, the Otosure is a PC-based automatic screening audiometer designed to provide occupational health professionals with significant workflow improvements through time and cost savings.

#### **Customised test options**

Automatic and manual test modes enable all patient hearing thresholds to be established and categorised.

The Otosure also has configurable options so the test can be customised to the users' specific requirements, ensuring accurate and efficient testing.

#### Data management

Offering seamless PC connectivity, the Otosure includes our Audibase data management software, allowing all test results to be stored and available for future review.

Extensive data and trend analysis provides business intelligence to support the delivery of a successful hearing conservation programme and identify at risk personnel.

As well as this, the instrument offers custom built educational tools with a predictive hearing level and hearing loss indicator, providing a clear explanation to test subjects, and enabling a greater level of assessment. All tests can be automatically categorised in accordance with HSE standards.

#### **EMR** connectivity

The Otosure is configured for use with the Audibase data management software and single-click test initiation.

The audiometer can also be interfaced to a number of leading OH Electronic Medical Record (EMR) systems to provide a paperless and seamless health screening programme that maximises workflow efficiencies.

#### **Portability**

Pocket-sized and weighing only 165g (0.36lbs), the Otosure is completely portable and comes supplied with a custom carry case, making it an ideal choice for multi-site audiometry testing.

#### Key features

Automatic and manual test modes

Automatic HSE categorisation

Unique educational and assessment tools

Extensive data analysis and trending

Third-party EMR connectivity

Compact, lightweight and portable

### Otosure

#### PC-based automatic screening audiometer

#### Air conduction audiometry

Air conduction range (kHz):	0.125, 0.25, 0.5, 1, 1.5, 2, 3, 4, 6, 8
Frequency accuracy:	<1%
Distortion:	<2%
Output level range:	-10dBHL to 100dBHL ±3dB
Output level step size:	Computer: 5dB Manual: 5dB
Test method:	Manual and automatic audiometry PC controlled, Computer (Hughson Westlake, BS6655)
Threshold level:	Computer: 5dB Manual: 5dB

#### Standard equipment

- Standard audiometric headset
- Patient response switch
- Carry case
- USB stick (manuals and Audibase software)
- USB cable (PC connection)

#### Optional equipment

- Amplivox Audiocups (noisereducing enclosures)
- ER-75 electro acoustic ear simulator
- Audiology booth
- Ear cushion covers for both standard and audiocup type headsets
- Booth leads


#### System requirements (PC operated)

Operating system:	Stand alone PC with Windows 10 or higher
Memory:	Internal memory and available disc space as required by the PC operating system
Software:	Supplied via USB stick
Interface:	USB

#### Physical data

Power:	Via USB connection from PC
Dimensions (L x W x H):	120 x 86 x 28mm
Weight:	165g / 0.36lbs

#### Additional information

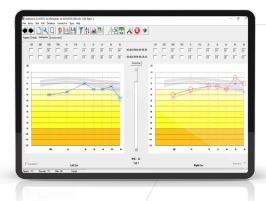

#### Safety and standards

Туре:	Audiometer Type 4
Safety:	IEC 60601-1 (plus UL, CSA & EN deviations)
EMC:	IEC 60601-1-2
Performance:	IEC 60645-1/ANSI S3.6
CE Mark:	Complies to EU Medical Device Regulation (MDR 2017/745)

#### **Audibase**

Developed by occupational health specialists for occupational health specialists, Audibase is designed to both manage and provide detailed information on hearing conservation programme performance.

Audibase can provide several features to assist with the interpretation of audiometric data, including automatic audiogram categorisation, user selected predictive displays (based upon age and gender) and hearing level indicators.



Amplivox Ltd, 3800 Parkside, Solihull Parkway, Birmingham Business Park, Birmingham, West Midlands, B37 7YG, United Kingdom

www.amplivox.com | +44 (0)1865 880846 | hello@amplivox.com